Soozier

IN230500580V02_US_CA

A90-290

EN

IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY

ASSEMBLY & INSTRUCTION MANUAL

Before using the product, please read the following instructions carefully. If you have the following conditions, please pay attention to your exercise time and consult a doctor before use.

- * Have heart disease
- * Have sensory impairment
- * Suffering from osteoporosis or diabetes
- * During pregnancy
- * Have a malignant tumor
- * Physical disability
- * Have high blood pressure or acute illness
- * Suffering from ventilation and arthritis
- * There are symptoms of dizziness, tinnitus, it is recommended to consult a doctor before using the product.

Precautions

- * Training against the wall

If the user can't control the direction of the force is downward, the training vehicle will move when used, and it needs to be used against the wall.


- * Seat height

Please choose a suitable chair for training, so as to avoid the phenomenon of unstable base caused by too high and too low

Safety warning

- 1.Do not stand on the pedal.
- 2.Make sure to sit on a chair, sofa or bed when exercising.
- 3.Before using the product, please make sure that all screws, bolts, knobs, frames, and brackets have been fixed.
- 4.Patients, the disabled and the elderly need to be used at the same time.
- 5.Those who suffer from heart disease or who are prohibited from doing strenuous exercise, please follow your doctor' s advice to exercise properly.
- 6.Please do not use this product for injured or unhealed
- 7.Please be sure to supervise children and the disabled when using it.
- 8.This product is designed for home use

PRODUCT ALL PARTS LIST

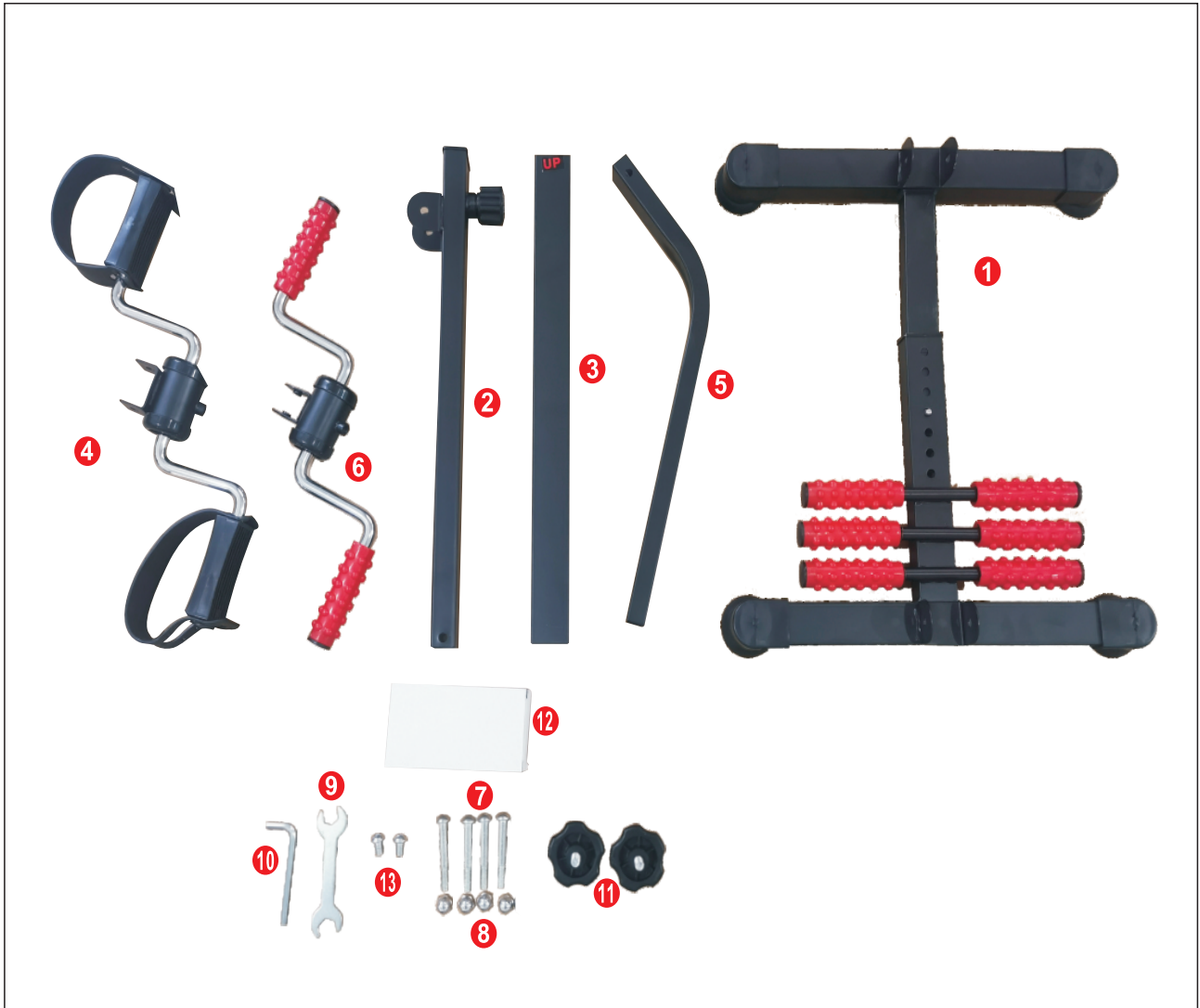
Parts list	Quantity
Lower main frame①	1
Central main frame②	1
Middle inclined connecting frame③	1
Lower limb rotation pedal assembly④	1
Upper connecting elbow⑤	1
Whole hand parts ⑥	1
Fixing screw⑦	4

Parts list	Quantity
Fixed nut ⑧	4
Multi-purpose wrench ⑨	1
Allen wrench ⑩	1
Plum knob ⑪	2
Meter ⑫	1
Fixing screw⑬	2

PRODUCT INSTALLATION STEPS

Step 1: First take out all the accessories in the carton and lay it down

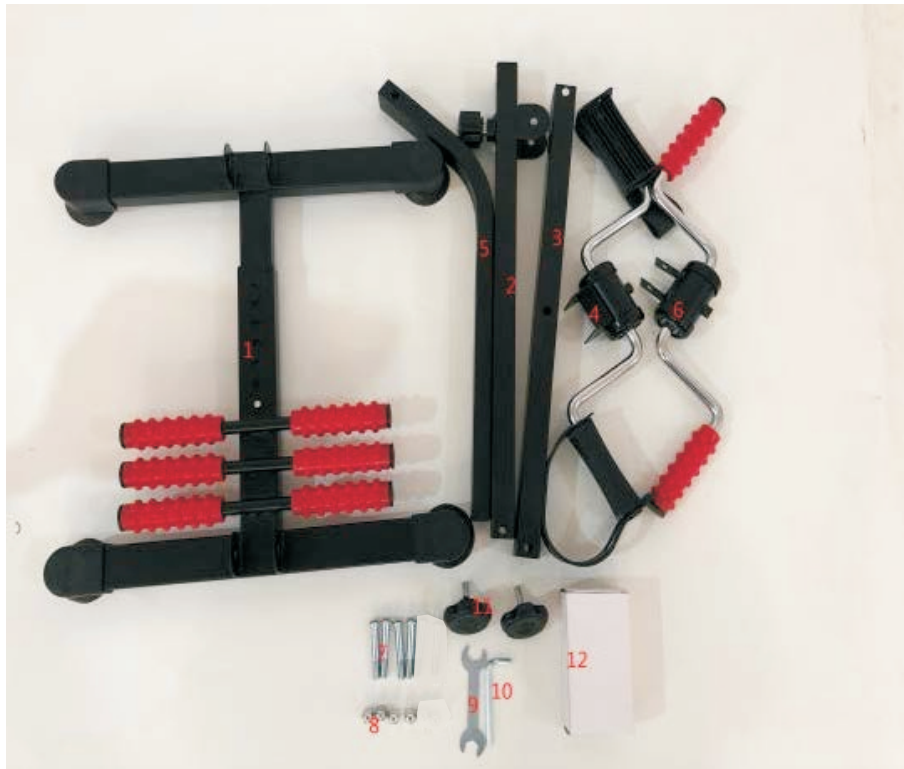Step 2: Lay the lower main frame ① flat and overlap the pipe hole on the left ②, ⑦ 8*50 screw and ⑧ use wrench ⑨ ⑩ to tighten together.

Step 3: Overlap the lower main frame ① with the right ③ pipe hole, and tighten the 8*50 screw and ⑧ with a wrench ⑨ ⑩.

Step 4: Overlap the fixed pipe holes ② and ③, and tighten the screws ⑦ 8*50 and ⑧ with wrench ⑨ ⑩.

Step 5: Splice the lower limb rotating foot assembly ④ on the upper end of the lower main frame ③, and insert the ⑦ 8*50 screw and ⑧ with the wrench ⑨ ⑩ to tighten the double holes.

Step 6: Insert the upper main frame into the upper end of the lower main frame, and then use the knob to pull the pin to fix it, adjust the upper and lower levels adaptively and flexibly, and then tighten the knob.

Step 7: Join the upper body rotating armrest assembly ⑥ at the upper end of the upper main frame , and insert the ⑬ screw with the wrench ⑨⑩ to tighten the double holes.

Step 8: The installation is complete!

Note: The electronic watch ⑫ can also be installed on the upper limb rotation pedal assembly ⑥ at your choice.

(18*50 screw is a spare screw)

US

If you have any questions, please contact our customer care center.
Our contact details are below:

001-877-644-9366

customerservice@aosom.com

Imported by Aosom LLC
27150 SW Kinsman Rd Wilsonville, OR 97070 USA
MADE IN CHINA

CA

If you have any questions, please contact our customer care center.
Our contact details are below:

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416-792-6088

customerservice@aosom.ca

Imported by Aosom Canada Inc.
7270 Woodbine Avenue, Unit 307, Markham, Ontario Canada L3R 4B9
MADE IN CHINA

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