
Soozier

IN240800132V01_US_CA

A90-385V00

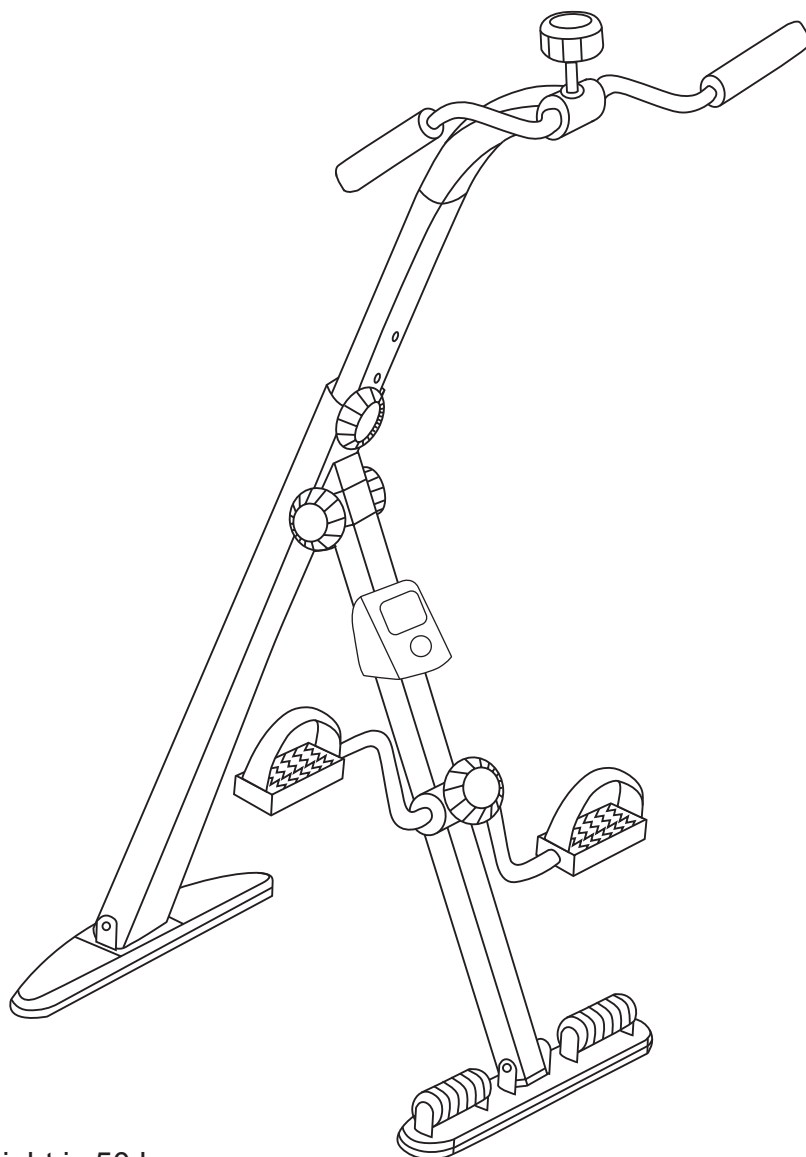* Maximum user weight is 50 kg.

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

EN: IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

① x 1

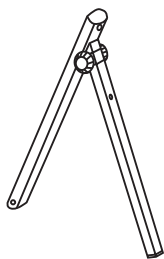
② x 1

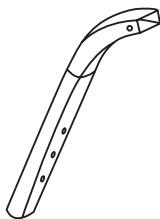
③ x 1


④ x 1

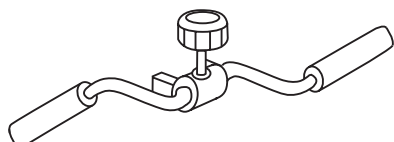
⑤ x 1

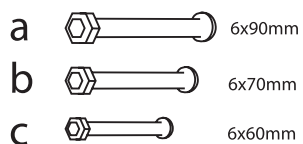
⑥ x 1

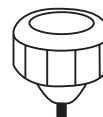
⑦ x 1

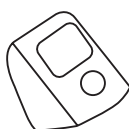
⑧ x 1

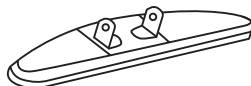
⑨ x 1

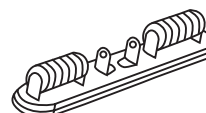⑩ x 2

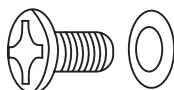⑪ x 1



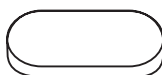
⑫ x 1

⑬ x 1

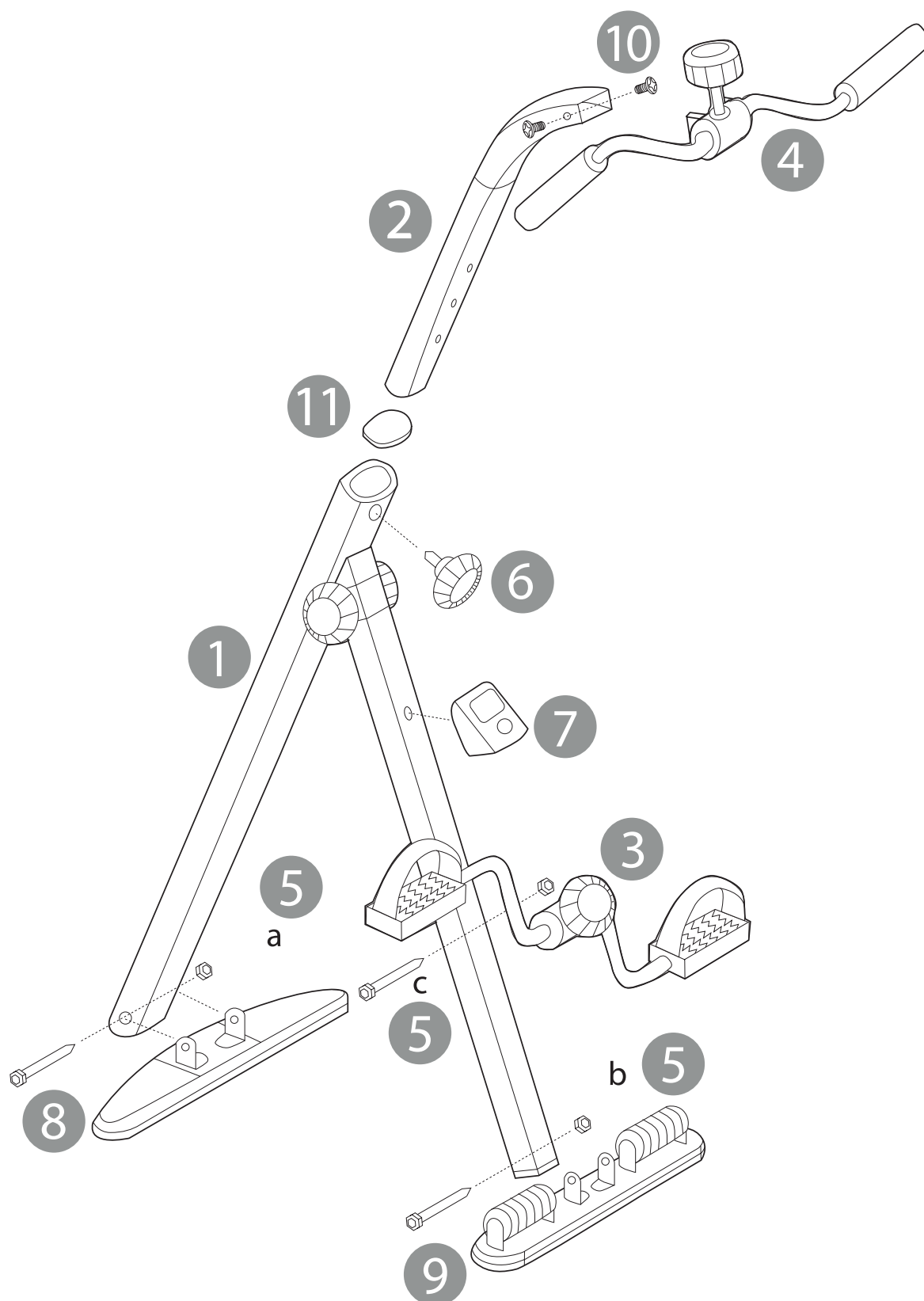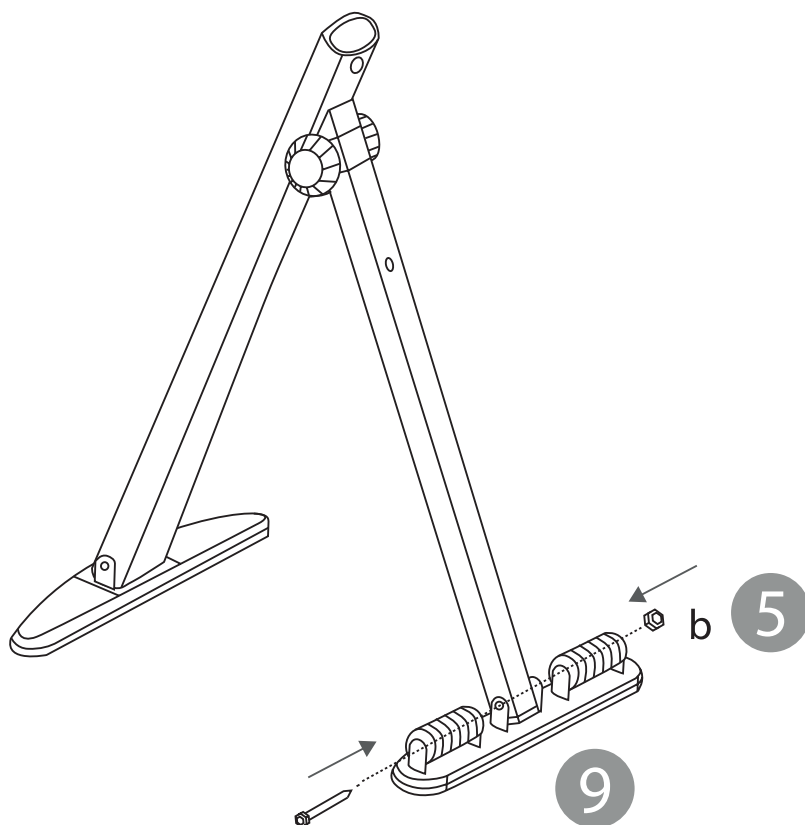
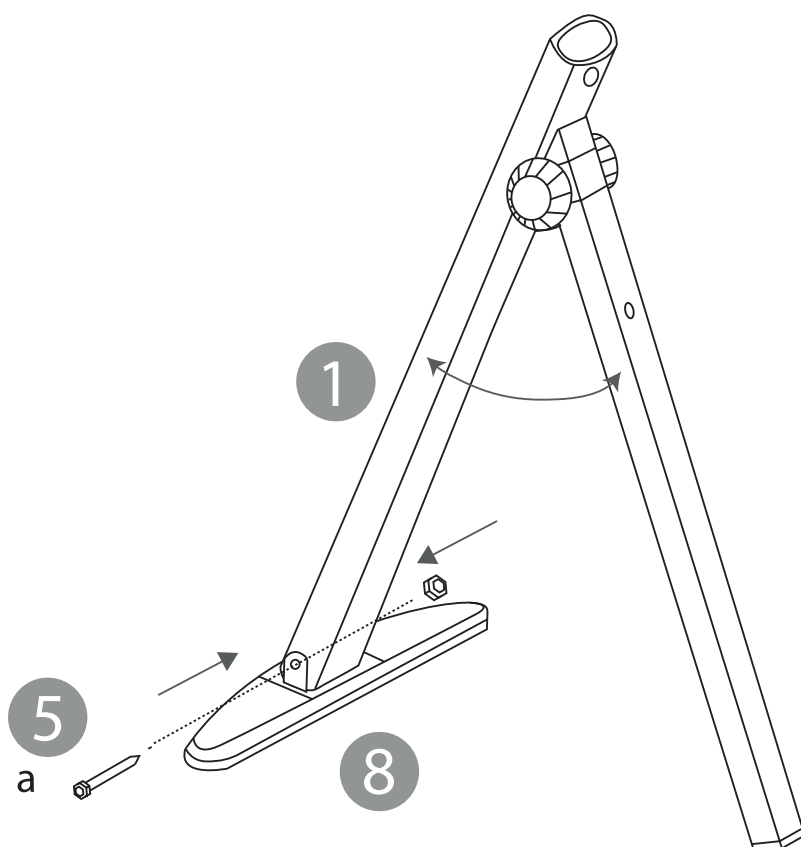
The magnet found in the pedals should be aligned to the right, with the rear pitchfork's sensor.

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
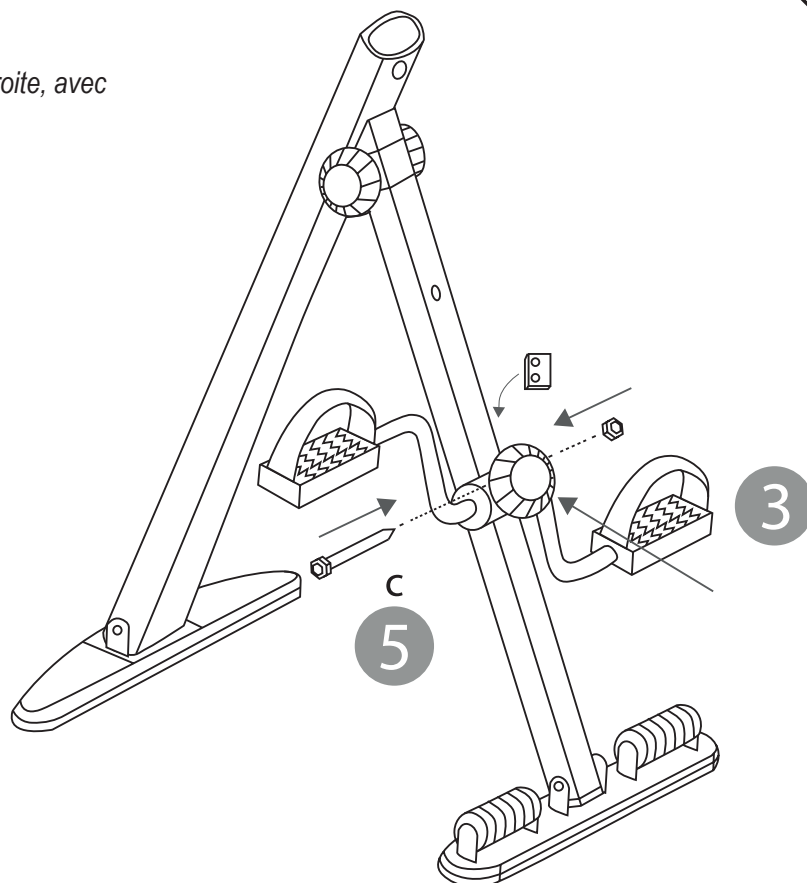
Insert the top front fork into the front fork and stop at the best position for you by tightening the adjustment knob (In case you only want to use the pedals without using the upper part, please place the plastic lid (11) on the hole for the lower pitchfork).

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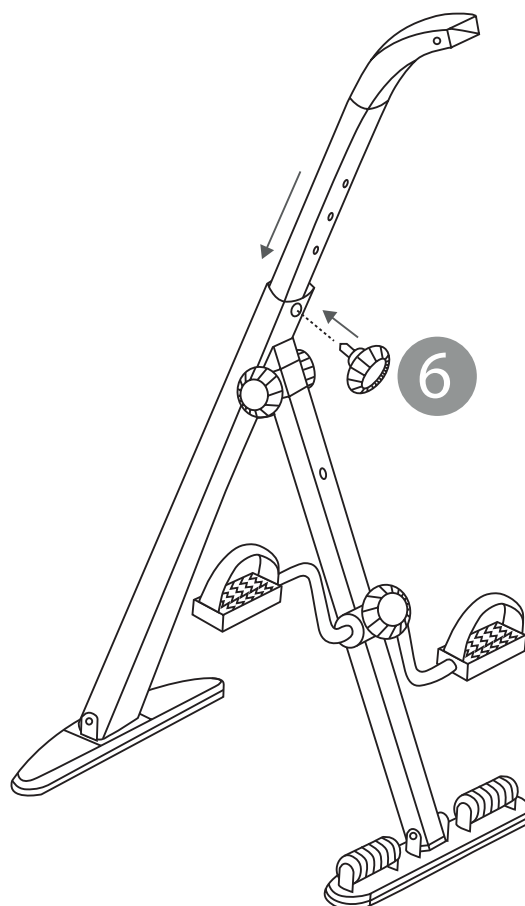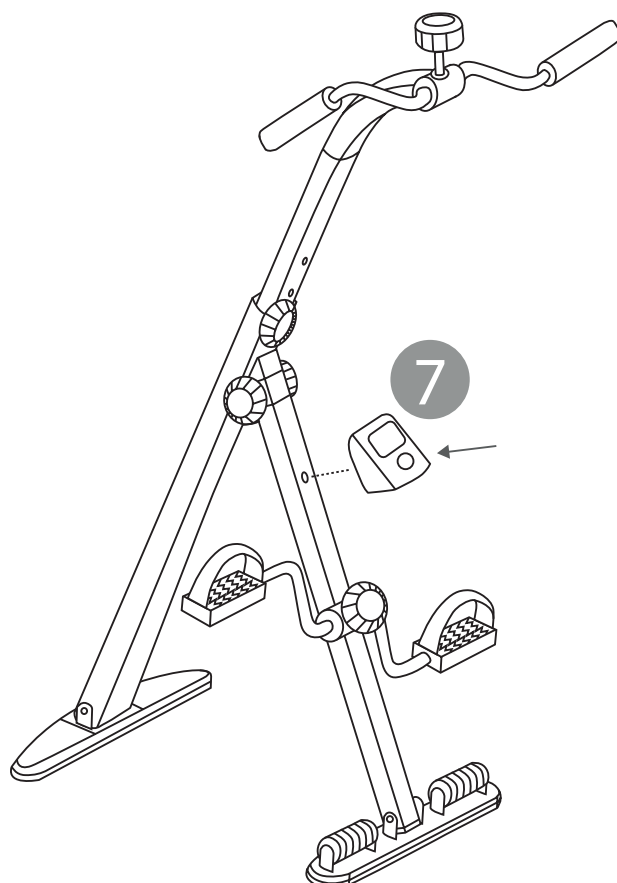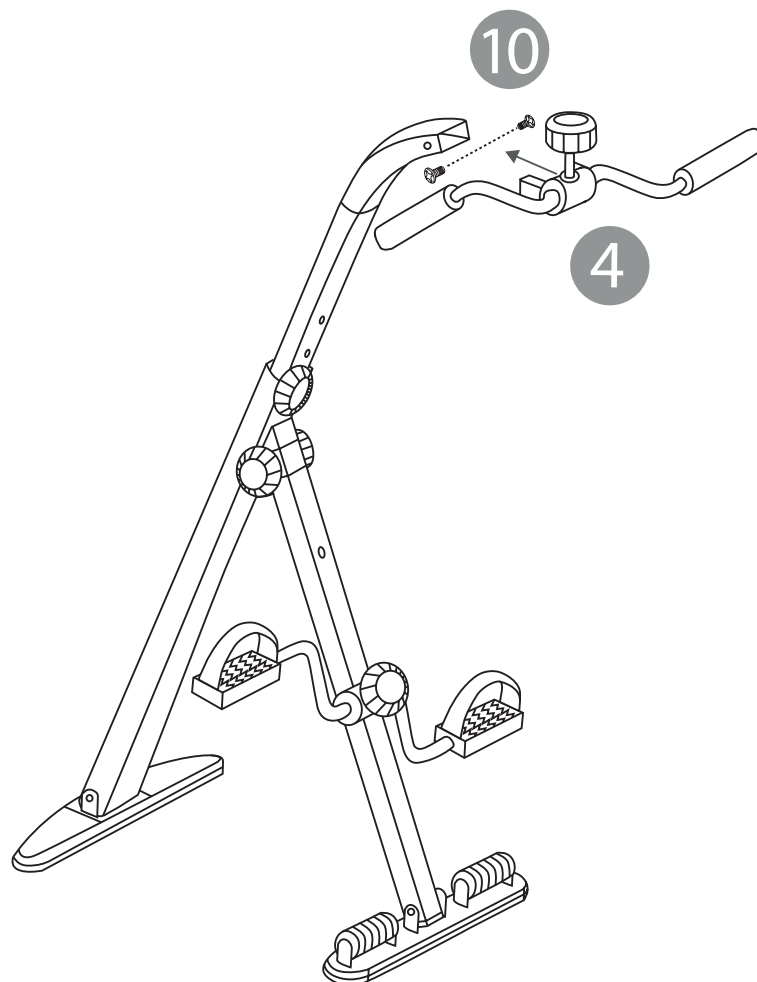READING CAREFULLY!**Please keep this manual in a safe place for easy reference**

1. It is important to read this entire manual before assembling and using the equipment. Safe and effective use can only be achieved if the equipment is assembled, maintained and used properly. It is your responsibility to ensure the all users of the equipment are informed of all warnings and precautions.
2. Before starting any exercise program you should consult your doctor to determine if you have any medical or physical conditions that could put your health and safety at risk, or prevent you from using the equipment properly. Your doctor's advice is essential if you are taking medication that affects your heart rate, blood pressure or cholesterol level.
3. Be aware of your body's signals. Incorrect or excessive exercise can damage your health. Stop exercising if you experience any of the following symptoms: pain, tightness in your chest, irregular heartbeat, extreme shortness of breath, lightheadedness, dizziness, or feelings of nausea. If you do experience any of these conditions, you should consult your doctor before continuing with your exercise program only.
4. Keep children and pets away from the equipment. The equipment is designed for adult use only.
5. Use the equipment on a solid, flat level surface with a protective cover on your floor or carpet. To ensure safety, the equipment should have at least 0.5 meters of free space all around it.
6. Before using the equipment, check that the nuts and bolts are securely tightened.
7. The safety of the equipment can only be maintained if it is regularly examined for damage and/or wear and tear.
8. Always use the equipment as indicated. If you find any defective components while assembling or checking the equipment, or if you hear any unusual noises coming from the equipment during use, stop immediately. Do not use the equipment until the problem has been rectified."
9. Wear suitable clothing while using the equipment. Avoid wearing loose clothing that may get caught in the equipment or that may restrict or prevent movement.
10. Suitable for home use only. Braking ability is independent of speed.
11. The equipment is not suitable for therapeutic use.
12. Care must be taken when lifting or moving the equipment so as not to injure your back. Always use proper lifting techniques and/or seek assistance necessary.

MAINTAIN

Neutral detergent can be used to dilute the water to wipe the cloth carefully wipe the oil, and then dry with a dry cloth.

Do not use alkaline solvents, solvents such as solvent, volatile oils and other solvents.

MAINTENANCE / INSPECTION

Avoid high temperature, moisture or direct sunlight, please place in a well-ventilated place.

When not in use please fold the admission, and placed in the children can not touch the Department.

If you do not use it for a while, be sure to check the precautions before using this product and use it for normal use.

When you accidentally wet the product, use a dry cloth to avoid skidding.


In case of malfunction, abnormality or any operational problem, please contact us at the company.

PRODUCT RECYCLING


Due to the different national recycling regulations, please refer to the relevant regulations.

INSTALLATION OF MONITOR


1. Insert the male wire plug from the rear fork into the female wire plug on the back of the monitor.

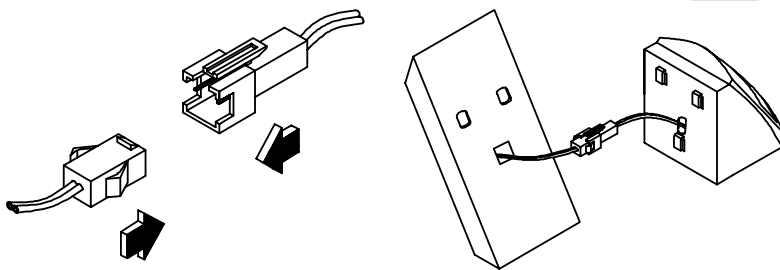2. Insert the joint wire into the rear fork.

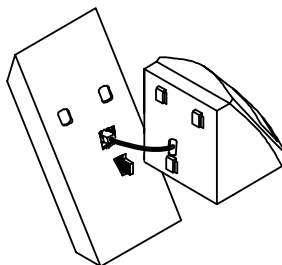
3. Slide the 3 hooks on the back of the monitor into the 3 outlets on the rear fork.

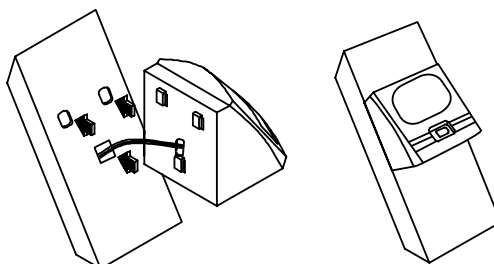
BATTERY REPLACEMENT


1. Open the top lid of the monitor using a small screwdriver.

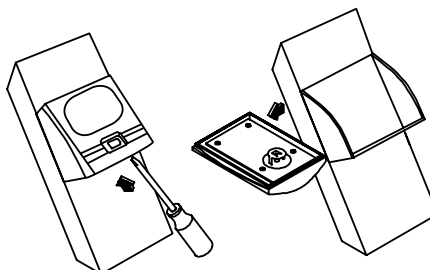
2. Remove the LR44 battery from its location.

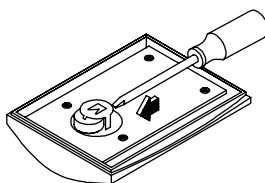3. Place the new LR44 battery [Positive (+) Side up] into position.

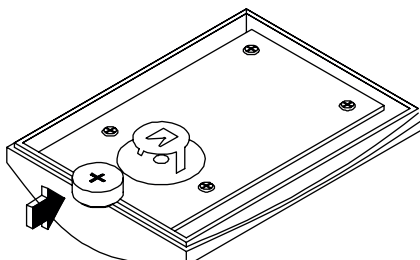4. Press the top lid of the monitor to close.

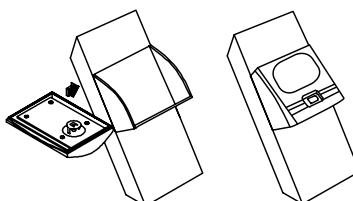USER MANUAL FOR THE MONITOR. TYPES OF FEATURES

CNT feature (number of pedalling revolutions): This feature allows us to count the number of cycling revolutions that we perform on every use of the product.

TM feature (time): This feature allows us to time how long we use our product. When we stop pedalling, the timer stops and continues the countdown when we start up once again.

CAL feature (calories): This feature allows us to count the number of calories that we burn by using the product.

DIST feature (distance): This feature allows us to know the distance travelled with the product. We travel 0.01 km for every 50 strokes we perform with the pedal.

TCNT feature (total pedalling counter): This feature allows us to count the number of strokes performed with the pedals starting from the first time we use the product. This function can only be reset by changing the monitor's battery.

STOP feature: The STOP signal appears when we stop pedalling and disappears when go back to putting our product in motion.

SCAN feature: By selecting the SCAN option, the CNT / TM/ CAL / DIST / TCNT features appear alternately on the monitor every 4 seconds.

To reset and put all the markers back to 0, press the red button on the monitor for 4 seconds. That will reset all the features except the TCNT feature, which will only reset by changing the monitor's battery.

HOW TO USE:

1. To turn on the monitor, push the red button. The counter at 0 and the STOP signal will appear.
2. Push the red button, as many times as necessary, until we enter into the feature that we want to see on the monitor and then start pedalling. While performing the exercise, we can push the red button to change the feature and see our progress or go to the SCAN mode so that the different features appear alternately and automatically on the monitor.
3. At the end of the exercise, the monitor automatically turns off after 4 minutes without using the product.

BATTERY INSTALLATION

The monitor operates with a coin cell battery. Under conditions of normal use, the battery can last up to six months. To change the battery, open the monitor cover, remove the used battery and insert a new one.

WARNING

- **INGESTION HAZARD:** This product contains a button cell or coin battery.
- **DEATH** or serious injury can occur if ingested.
- A swallowed button cell or coin battery can cause **Internal Chemical Burns** in as little as **2** hours.
- **KEEP** new and used batteries **OUT OF REACH of CHILDREN**
- **Seek immediate medical attention** if a battery is suspected to be swallowed or inserted inside any part of the body.

1.) Battery type: LR44

2.) Nominal Voltage 1.5V

Battery

- Remove and immediately recycle or dispose of used batteries according to local regulations and keep away from children. Do NOT dispose of batteries in household trash or incinerate.
- Even used batteries may cause severe injury or death.
- Call a local poison control center for treatment information.
- A statement indicating the compatible battery type LR44.
- A statement indicating the nominal battery voltage.
- Non-rechargeable batteries are not to be recharged.
- Do not force discharge, recharge, disassemble, heat above (manufacturer's specified temperature rating) or incinerate. Doing so may result in injury due to venting, leakage or explosion resulting in chemical burns.
- Ensure the batteries are installed correctly according to polarity (+ and -).
- The statement "Do not mix old and new batteries, different brands or types of batteries, such as alkaline, carbon-zinc, or rechargeable batteries."
- The statement "Remove and immediately recycle or dispose of batteries from equipment not used for an extended period of time according to local regulations."
- The statement "Always completely secure the battery compartment. If the battery compartment does not close securely, stop using the product, remove the batteries, and keep them away from children."

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Imported by Aosom LLC
Imported by Aosom Canada Inc.
Importé par Aosom Canada Inc.